

1 **Article type:** Original Research Article

2

3 **Running headline:** Highly sensitive virus detection from saliva

4

5 **Title:** Development of a simple and highly sensitive virion concentration method to
6 detect SARS-CoV-2 in saliva

7

8 **Authors:** Yasuko Yamazaki, Uxía Alonso Alonso, Remil L. Galay, Wataru Yamazaki

9

10 **Affiliations:**

11 Center for Southeast Asian Studies, Kyoto University, 46 Shimoadachi-cho, Yoshida,

12 Sakyo-ku, Kyoto 606-8501, Japan (Y. Yamazaki, yasuko.y@cseas.kyoto-u.ac.jp; U.

13 Alonso Alonso, alonsouxia93@gmail.com; R. L. Galay, rlgalay@up.edu.ph; W.

14 Yamazaki, yamazaki@cseas.kyoto-u.ac.jp)

15 Animal Medicine and Health at Institut de Recerca i Tecnologia Agroalimentàries

16 (IRTA) – Centre de Recerca en Sanitat Animal (CRESA), Spain (U. Alonso Alonso)

17 Department of Veterinary Paraclinical Sciences, College of Veterinary Medicine,

18 University of the Philippines Los Baños, Los Baños, Laguna 4031, Philippines (R. L.

19 Galay)

20 Kyoto University School of Public Health, Konoe-cho, Yoshida, Sakyo-ku, Kyoto 606-

21 8303, Japan (W. Yamazaki).

22

23 **Address for correspondence:**

24 Wataru Yamazaki, Center for Southeast Asian Studies, Kyoto University, 46

25 Shimoadachi-cho, Yoshida, Sakyo-ku, Kyoto 606-8501, Japan; email:

26 yamazaki@cseas.kyoto-u.ac.jp

27

28 **Abbreviations:**

29 ATCC, American Type Culture Collection

30 LOD, Limit of detection

31 PBS, Phosphate-buffered saline

32 PEG, Polyethylene glycol

33 RT-qPCR, Reverse transcription quantitative real-time polymerase chain reaction

34 SAP, Semi-alkaline proteinase

35 SARS-CoV-2, Severe acute respiratory syndrome coronavirus 2

36

37 **Highlight:**

38 ☐ A method has been developed to detect SARS-CoV-2 from human saliva with
39 100 times higher sensitivity than conventional methods.

40 ☐ The developed method combines simple pretreatment within 60 min with
41 conventional nucleic acid extraction and RT-qPCR.

42 ☐ This method can be applied for more sensitive virus testing from individual saliva.

43 ☐ This method can potentially be applied to screening more than 100 saliva
44 samples while maintaining the equivalent detection power of conventional
45 methods.

46 ☐ The method can be adapted to improve the sensitivity of detecting various
47 pathogens from human and animal saliva.

48

49 **Abstract**

50 **Background:** Controlling novel coronavirus pandemic infection (COVID-19) is a
51 global challenge, and highly sensitive testing is essential for effective control. The
52 saliva is a promising sample for high-sensitivity testing because it is easier to collect
53 than nasopharyngeal swab samples and allows large-volume testing.

54 **Results:** We developed a simple SARS-CoV-2 concentration method from saliva
55 samples that can be completed in less than 60 min. We performed a spike test using
56 12 ml of saliva samples obtained from healthy volunteer people, and the developed
57 method performance was evaluated by comparison using a combination of automatic
58 nucleic acid extraction followed by RT-qPCR detection. In saliva spike tests using a
59 10-fold dilution series of SARS-CoV-2, the developed method was consistently 100-
60 fold more sensitive than the conventional method.

61 **Conclusions:** The developed method can improve the sensitivity of the SARS-CoV-
62 2 test using saliva and speed up and save labor in screening tests by pooling many
63 samples. Furthermore, the developed method has the potential to contribute to the
64 highly sensitive detection of various human and animal viral pathogens from the
65 saliva and various clinical samples.

66

67 **Keywords:** Concentration, COVID-19; SARS-CoV-2, Semi Alkaline Proteinase;
68 virion, virus

69

70 Introduction

71 The COVID-19 epidemic continues as of 2023 and remains a public health threat
 72 (WHO). Virus detection using purified RNA obtained by extraction kits and real-time
 73 reverse transcription quantitative PCR (RT-qPCR) has been universally used due to
 74 its high detection sensitivity and low incidence of false negatives due to nonspecific
 75 amplification (Lu *et al.*, 2020; Vogels *et al.*, 2020). Virus detection using saliva has
 76 been used to diagnose respiratory infections such as COVID-19 and influenza
 77 because of its easy sampling and low burden on patients (Azzi *et al.*, 2020; To *et al.*,
 78 2017; Vogels *et al.*, 2020). However, because of its somewhat lower detection
 79 sensitivity compared to nasopharyngeal swab samples, false-negative results occur
 80 in samples collected early in infection or late in recovery when viral load is low,
 81 meaning that patients who slip through the test may not receive appropriate
 82 quarantine measures and become a potential source of infection (Azzi *et al.*, 2020;
 83 To *et al.*, 2017; Vogels *et al.*, 2020; Yamazaki *et al.*, 2021).

84

85 In addition, in the early stages of an epidemic, when there are few positive
 86 patients, it is crucial to test a large number of samples from a large number of people
 87 for negative confirmation to prevent the spread of infection (Barat *et al.*, 2021;
 88 Watkins *et al.*, 2021). In this case, to save cost and labor, pooling of samples, such
 89 as saliva (Barat *et al.*, 2021; Watkins *et al.*, 2021), and nasal, nasopharyngeal and
 90 oropharyngeal swabs (Ayaz *et al.*, 2022; Praharaj *et al.*, 2020; Pratelli *et al.*, 2022)
 91 from several people and testing them together is sometimes used for screening.

92

93 However, since the positive samples could be diluted by mixing with negative
 94 samples, leading to low virus concentration in the pooled sample, the test can

become a false negative if it is below the detection limit; hence, accurate detection may not be possible. As a solution, we have developed a method for detecting concentrated viruses in samples through immunomagnetic beads (Yamazaki *et al.*, 2019; Makino *et al.*, 2020). Still, it is not versatile because it requires specific antibodies for each virus.

The Polyethylene glycol (PEG) precipitation method has been used worldwide to enrich and detect norovirus and other viruses from oysters (Lowther *et al.*, 2019; National Institute of Health Sciences, Japan (NIHS) 2010; Yamazaki *et al.*, 2022). While this method can concentrate any virus, the presence of sample-derived inhibitors reduces concentration performance (Lowther *et al.*, 2019; Miura *et al.*, 2018; Yamazaki *et al.*, 2022). In our previous studies, we have shown that a combination of a very short, low centrifugation process (900 *g*, 1 min) and a normal centrifugation process (10,000-20,000 *g*, 5 min) as a pretreatment step for simple concentration detection of target bacteria in chicken cecal contents (Sabike *et al.*, 2016). Also, we have demonstrated that genetic testing for SARS-CoV-2 is possible without using an extraction kit by digesting human saliva containing potential genetic testing inhibitors with semi-alkaline protease (SAP) (Yamazaki *et al.*, 2021). Here, we report the successful development of a new method for the concentration and detection of SARS-CoV-2 from a large volume of saliva by improving and integrating our previously published methods.

Materials and Methods

Saliva sampling

Saliva samples were collected from three healthy volunteers, i.e., the three authors

of this paper (YY, UAA, and WY), by repeatedly transferring drool collected in the oral cavity into a 50-ml sterile tube. After each saliva was thoroughly mixed by vortexing, the multiple saliva was promptly mixed in a new 50-ml tube to produce approximately 49 ml of pooled saliva. The three saliva samples used were confirmed to be SARS-CoV-2 negative by two RNA extraction methods (conventional and developed) and RT-qPCR detection, as described below, before the experiment.

Preparation of a 10-fold dilution series of SARS-CoV-2 spiked saliva

A 10-fold dilution series of heat-inactivated SARS-CoV-2 (ATCC VR-1986HK; American Type Culture Collection, Manassas, VA, USA) in PBS was prepared. The pooled saliva was dispensed into four 50-ml tubes of 12.2 ml each. SARS-CoV-2-containing saliva from neat to 10(-4) fold dilutions was prepared by sequentially spiking the 10-fold dilution series of SARS-CoV-2 into the 50 ml tubes containing pooled saliva and then vortexed thoroughly (Table 1).

RNA extraction by the conventional method

According to the pathogen detection manual 2019-nCoV issued by the National Institute of Infectious Diseases, Japan (NIID-J), two sets of the 100 µl saliva-spiked SARS-CoV-2 were collected into 1.5-ml microcentrifuge tubes and diluted 1:3 with 300 µl of PBS and sputum homogenizer SAP (Semi-Alkaline Proteinase, Suputazyme; Kyokuto Pharmaceutical Industrial, Tokyo, Japan), respectively. After sufficient vortexing, both (PBS and SAP) of the 1:3 dilutions containing SARS-CoV-2 were centrifuged at 20,000 *g* for 30 min, the former immediately and the latter after 15-min incubation with ten manual inversions mixing every 3 min at room temperature. The resulting 200 µl of the supernatant was transferred in a new 1.5-ml

microcentrifuge tube and was set in an automated nucleic acid extractor MagLead 6GC (Precision System Science, Co., Ltd, Matsudo, Japan) with MagDEA Dx SV reagent cartridge (Precision System Science) and were extracted and purified as RNA in 50 µl of distilled water.

Viral concentration and RNA extraction by developed method

An overview is shown in Figure 1. Specifically, 12 ml of SAP (Kyokuto) was added to the remaining 12 ml of saliva containing SARS-CoV-2. After vortexing, the mixture was kept at room temperature for 15 min. During the 15-min incubation, ten inversion mixings were performed manually every 3 min. Then, 4,000 *g*, 5 min initial centrifugation was performed. Taking care not to inhale the pellet derived from the formed saliva components, 18 ml (75% of the initial mixture volume) of the centrifugal supernatant was prudently transferred to a new 50-ml tube. 15 ml of SAP (Kyokuto) was added and incubated for 15 min at room temperature. Then, 13.2 ml of PEG solution (40% PEG-NaCl, see details in our previous publication, Yamazaki *et al.*, 2022) was added and thoroughly mixed by vortexing, followed immediately by a second centrifugation at 8,000 *g* for 20 min. The supernatant was carefully removed after the second centrifugation. To prevent contamination, 100 µl of PBS in a 1-ml long tip, which is longer than the 50 ml tube, was added. Pipetting was performed with the 1-ml long tip set pipet from the bottom of the tube to the area where the pellet had adhered during the first centrifugation, where precipitates of PEG-virus particle complex are assumed to be attached, although it cannot be seen with the naked eye. In addition, to completely detach any PEG-virus particle complex precipitates that may still be adhering to the tube wall, the 1-ml short tip was added to the 50-ml tube, the lid was closed, and the tube was thoroughly vortexed.

Approximately 200 µl of the mixture of about 100 µl of PBS (containing PEG-virus particles) added to the around 100 µl of supernatant remaining on the inner wall of the 50-ml tube obtained by flushing was transferred into a 1.5-ml screw cap tube using a 1-ml long tip. The mixture was then extracted and purified as 50 µl of RNA using an automated nucleic acid extractor (Precision System Science), as described above. When the extracted RNA could not be tested immediately, it was stored at -80°C until use.

Conducting RT-qPCR and determination of LOD

RT-qPCR was performed with 4 µl of the extracted RNA in 20 µl of the reaction mixture using a QuantStudio 3 (Thermo Fisher Scientific, Inc., Waltham, MA, U.S.A.), according to the method by Lu and colleagues (2020). The amplification time was slightly extended to ensure detection, as described below. Details of the reagents used are as follows: 20-µl RT-qPCR reactions comprised 10 µl of SuperScript III Platinum One-step RT-qPCR 2x reaction (Thermo Fisher Scientific), 0.4 µl of SuperScript III/Platinum *Taq* Mix (Thermo Fisher Scientific), 2 µl of primer (Hokkaido System Science Co. Ltd., Sapporo, Japan) probe (Integrated DNA Technologies, Inc, Singapore) mix for SARS-CoV-2 N2 detection reported by Emergency Use Authorization issued by the US Food and Drug Administration issued by the FDA (Lu *et al.*, 2020), 2 µl of magnesium sulfate (50 mM, Thermo Fisher Scientific), 0.2 µl of Rox Dye (Thermo Fisher Scientific) diluted 1:5 with distilled water, 1.4 µl of nuclease-free water, and 4 µl of the RNA template. The cycling conditions were as follows: one cycle at 50°C for 900 sec and 95°C for 120 sec, followed by 50 cycles each at 95°C for 15 sec and 55°C for 60 sec. The automatically calculated Ct value was adopted, and the Ct cut-off value was set at

40.000. Positive results were determined if all three Ct values were within 40.000 in triplicate. Samples that showed only one or two positive Ct values in the triplicate analysis were interpreted as negative. The mean and standard deviation of the Ct were calculated for all samples interpreted as positive.

Results

Tables 1 and 2 show that the developed concentration method enabled 100-fold more sensitive detection than the conventional method by adding only a simple pretreatment within 60 min before the conventional extraction method. The developed method required only a centrifuge machine for 50 ml tubes up to 8,000 *g* and a vortex for mixing the liquid in the 50 ml tubes and did not need expensive equipment such as an ultracentrifuge.

As shown in Tables 1 and 2, the conventional method required 20,000 to 40,000 copies of virus per ml of saliva to detect SARS-CoV-2, whereas the developed method required only 200 copies of virus per ml of saliva. In other words, the developed method was at least 100 times more sensitive than the conventional method. Furthermore, a comparison of the number of viral copies per RT-qPCR reaction tube showed that 32 to 64 and 58 copies were required for detection by the conventional and developed methods, respectively. Namely, LOD per reaction tube was comparable for the two methods.

Discussion

In Japan, the airport quarantine for COVID-19 recommends collecting and submitting approximately 5 ml of saliva. Still, following the protocol of the NIID-J

(2020), 200 µl of the supernatant is generally used after diluting 200 µl within the range of 1:1 to 1:3 ratios and centrifuging at 20,000 *g* for 30 min. In other words, only 50-100 µl of saliva is used for the testing, and the remaining saliva of around 5 ml is used for nothing but retests. If all the unused samples were simply submitted to the concentration method developed in this study, the detection sensitivity would be dramatically increased up to 100 times. Hence, a more accurate quarantine control measure would be possible.

The disadvantage of pool testing is decreased LOD (Barat *et al.*, 2021; Praharaj *et al.*, 2020; To *et al.*, 2017; Watkins *et al.*, 2021), but the developed concentration method can solve this problem. Hence, the developed method is ideal for labor-saving large-scale screening in the early stages of an outbreak when the positivity rate is low. When a new variant emerges in the future, the concentration method could be used as a large-scale screening test to reasonably enhance quarantine control measures and contribute to efficient epidemic control. The developed method is theoretically capable of 180-fold virus concentration (Figures 1 and 2). Since the actual measured value is about 100-fold (Tables 1 and 2), it should be noted that the recovery rate could be reduced by approximately 50%. In the conventional method, the 12-ml saliva sample is equivalent to 240 pooled samples of 50 µl of saliva per individual. Still, considering the recovery rate, a pool of more than 100 samples can be expected to have a detection sensitivity equivalent to or better than that of the conventional method.

In our previous study, an immunomagnetic bead method using specific antibodies was successfully used to detect influenza A viruses added to PBS, duck

feces, and chicken meat at 10- to 1,000-fold sensitive concentration (Yamazaki *et al.*, 2019; Makino *et al.*, 2020). Although this method is extremely sensitive, it is not very versatile because it requires the preparation of specific antibodies for each virus species and has the disadvantage that the LOD is reduced when the samples contain many inhibitory substances, such as components of the duck feces and chicken meat. In the present study, this problem has been successfully overcome by improving the pretreatment method with the combination of SAP, a sputum dissolving agent, and the PEG precipitation method to achieve highly sensitive SARS-CoV-2 concentration and detection in saliva.

The PEG precipitation method customarily includes an overnight 4°C incubation process for virion capture (Lowther *et al.*, 2017; Yamazaki *et al.*, 2022). However, the NIH protocol (2010) states that the virion-PEG complex can be recovered from the oyster midgut gland immediately by centrifugation at 8,000 *g* for 20 min without an overnight incubation. In the present study, we referred to this finding and confirmed that SARS-CoV-2 in saliva could be concentrated immediately after being centrifuged at 8,000 *g* for 20 min without needing overnight incubation before the centrifuge, as expected. This allowed us to establish a rapid protocol successfully. This study's limitations include the inability to evaluate SARS-CoV-2-positive clinical samples and the fact that, although the concentration process is simple, the number of steps involved requires care by the examiner to avoid contamination and laboratory infection.

When the PEG-NaCl solution is mixed with a liquid sample containing trace amounts of virions, PEG, a polymer, adsorbs water molecules, causing the virions to

aggregate. By centrifuging, the agglomerated virions-PEG complex can be precipitated on the wall of the tube. On the other hand, if the liquid sample contains impurities, this reaction is inhibited, and the recovery rate of virus particles is reduced. As shown in Table 1, the developed method is more sensitive than the conventional method, but the estimated number of copies of virus per reaction tube may be higher for the developed method to obtain a positive result. For example, as shown in Tables 1 and 2, Test 2 showed Ct values at 36.104 and 36.310 for the conventional method (64 copies/RT-qPCR reaction tube) versus the Ct value at 38.320 for developed method (115 copies/RT-qPCR reaction tube). This suggests that, although carefully collected, there is some residual material in the supernatant after centrifugation that inhibits virion-PEG complex formation or that there reduced viral recovery, resulting in an increase of the Ct value in the developed method.

In the developed method, only 75% of the saliva centrifugal supernatant is used and the remaining 25% must be discarded without use, as shown in Figure 1, Process 2. In our preliminary experiment, the maximum amount of centrifugal supernatant corresponding to more than 90% was collected, and concentration detection was attempted. However, contrary to expectations, the recovery rate was more than 10 times lower than the theoretical value (data not shown). In this case, a large pellet was identified on the tube by the naked eye in process 6 of Figure 1. In other words, although the centrifugal supernatant appeared clear by the naked eye observation in Processes 2-5, we speculated that saliva components were mixed in and inhibited the formation of virion-PEG complexes. Therefore, we did not use the lower portion of the centrifugal supernatant, which was presumed to contain more saliva components due to the gradient caused by centrifugation but used only the

upper 75% to obtain a stable recovery rate.

Nevertheless, the developed method has the potential to solve the technical limitation of conventional genetic testing methods, i.e., the problem that samples carrying trace amounts of the virus have been judged as false negative because they are below the LOD (Barat *et al.*, 2021; Praharaj *et al.*, 2020; To *et al.*, 2017; Watkins *et al.*, 2021; Yamazaki *et al.* 2019). For example, rabies is transmitted from dogs to humans via dog saliva. Still, definitive diagnosis requires dog brain emulsion containing large amounts of rabies virus, not the dog saliva, due to lack of the virus amount, and there is an animal welfare issue of euthanasia of dogs for sampling. The oral pulse oximeter used for anesthesia monitoring in veterinary clinics caused the transmission of severe fever thrombocytopenia syndrome (SFTS) virus transmission to cats due to contamination. However, the low amount of the virus was not detectable (Mekata *et al.*, 2023). Both cases demonstrate the problem that current genetic testing cannot accurately detect trace amounts of virus in saliva, leading to false-negative diagnoses.

Conclusions

We have successfully developed a simple and highly sensitive method for concentrating SARS-CoV-2 in saliva. We demonstrated that the developed method has an extremely sensitive detection performance that is at least 100 times higher than the conventional method. We further showed that this method has the potential to screen more than 100 saliva samples with power comparable to conventional extraction methods. In the future, this method may be applied to highly sensitive diagnosis of various human and animal viral infections that can be tested from saliva,

as well as to rapid screening by pooled testing of a large number of samples.

Ethics approval

The Ethics Committee of Kyoto University Graduate School and the Faculty of Medicine approved this study (R2379), and consent to participate was waived from three volunteer saliva sample donors who are also authors.

Declaration of Availability of data and materials

All data obtained in this study is included in the paper and Supplemental Table 1. In addition, the data sets in this study are available from the corresponding author upon reasonable request.

Declaration of Competing Interests

The authors declare that they have no known competing interests.

Declaration of Consent for publication

Consent for publication has waived the need to obtain informed consent from three volunteer saliva sample donors who are also authors.

Declaration of Generative AI and AI-assisted technologies in the writing process

During the preparation of this work, the authors used DeepL (DeepL SE, Cologne, Germany) to improve readability and language. After using this tool, the authors reviewed and edited the content as needed and took full responsibility for the content of the publication.

345

346 **Declaration of Submission**

347 The authors confirm that this manuscript or data has not been previously published
348 and is not being considered for publication elsewhere. The authors further confirm
349 that all authors have contributed to the study and have approved the final version.

350

351 **Funding**

352 This research was supported by AMED Number JP20he0622031, JSPS KAKENHI
353 Numbers JP21H03180, JP22K05950, JP22KK0097, JSPS Bilateral Program
354 Number JPJSBP120199944, and the Joint Usage/Research Center for Global
355 Collaborative Research, Center for Southeast Asian Studies, Kyoto University.

356

357 **Authors' contributions**

358 **YY:** Data curation, Formal analysis, Funding acquisition, Methodology, Resources,
359 Validation, Visualization, Writing - original draft, Writing - review & editing. **UAA:**
360 Methodology, Resources. **RLG:** Formal analysis, Writing - review & editing. **WY:**
361 Conceptualization, Data curation, Funding acquisition, Methodology, Project
362 administration, Resources, Supervision, Validation, Visualization, Writing - original
363 draft, Writing - review & editing.

364

365 **Reference**

- 366 1. Ayaz A, Demir AGO, Ozturk G, Kocak M. A pooled RT-PCR testing strategy for
367 more efficient COVID-19 pandemic management. Int J Infect Dis. 2022;116:1-6.
368 doi:10.1016/j.ijid.2021.12.328
- 369 2. Azzi L, Carcano G, Gianfagna F, Grossi P, Gasperina DD, Genoni A, et al.
370 Saliva is a reliable tool to detect SARS-CoV-2. J Infect. 2020;81:e45-e50.

- 371 <https://doi.org/10.1016/j.jinf.2020.04.005>
- 372 3. Barat B, Das S, De Giorgi V, Henderson DK, Kopka S, Lau AF, Miller T, Moriarty
373 T, Palmore TN, Sawney S, Spalding C, Tanjutco P, Wortmann G, Zelazny AM,
374 Frank KM. Pooled Saliva Specimens for SARS-CoV-2 Testing. J Clin Microbiol.
375 2021 Feb 18;59(3):e02486-20. doi: 10.1128/JCM.02486-20.
- 376 4. Lowther JA, Bosch A, Butot S, Ollivier J, Mäde D, Rutjes SA, Hardouin G,
377 Lombard B, In't Veld P, Leclercq A. Validation of EN ISO method 15216 - Part 1 -
378 Quantification of hepatitis A virus and norovirus in food matrices. Int J Food
379 Microbiol. 2019 Jan 2;288:82-90. doi: 10.1016/j.ijfoodmicro.2017.11.014.
- 380 5. Lu X, Wang L, Sakthivel SK, Whitaker B, Murray J, Kamili S, et al. US CDC
381 Real-Time Reverse Transcription PCR Panel for Detection of Severe Acute
382 Respiratory Syndrome Coronavirus 2. Emerg Infect Dis. 2020;26(8):1654-1665.
383 <https://doi.org/10.3201/eid2608.201246>
- 384 6. Mekata H, Umeki K, Yamada K, Umekita K, Okabayashi T. Nosocomial Severe
385 Fever with Thrombocytopenia Syndrome in Companion Animals, Japan, 2022.
386 Emerg Infect Dis. 2023 Mar;29(3):614-617. doi: 10.3201/eid2903.220720.
- 387 7. Miura T, Schaeffer J, Le Saux JC, Le Mehaute P, Le Guyader FS. Virus Type-
388 Specific Removal in a Full-Scale Membrane Bioreactor Treatment Process. Food
389 Environ Virol. 2018 Jun;10(2):176-186. doi: 10.1007/s12560-017-9330-4.
- 390 8. National Institute of Health Sciences (NIHS), Committee for Standard Test
391 Methods for Viruses in Foods. Concentration of viruses from bivalve mollusks
392 (oysters). August 20, 2010 (in Japanese).
393 http://www.nihs.go.jp/fhm/csvdf/kentest/csvdf001_wg_100820.pdf (accessed
394 December 15, 2023).
- 395 9. National Institute of Infectious Diseases, Japan (NIID-J), Pathogen detection

manual 2019-nCoV. March 19, 2020 (in Japanese).

<https://www.niid.go.jp/niid/images/lab-manual/2019-nCoV20200319.pdf>

(accessed December 15, 2023).

10. Praharaj I, Jain A, Singh M, et al. Pooled testing for COVID-19 diagnosis by real-time RT-PCR: A multi-site comparative evaluation of 5- & 10-sample pooling.

Indian J Med Res. 2020;152(1 & 2):88-94. doi:10.4103/ijmr.IJMR_2304_20

11. Pratelli A, Pellegrini F, Ceci L, et al. Severe acute respiratory syndrome

coronavirus 2 detection by real time polymerase chain reaction using pooling

strategy of nasal samples. Front Microbiol. 2022;13:957957. Published 2022 Jul

22. doi:10.3389/fmicb.2022.957957

12. Sabike, II, Uemura, R, Kirino, Y, Mekata, H, Sekiguchi, S, Farid, AS, Goto, Y,

Horii, Y, Yamazaki, W. Assessment of the *Campylobacter jejuni* and *C. coli* in

broiler chicken ceca by conventional culture and loop-mediated isothermal

amplification method. Food Control, 74, 2017,107-111.

<https://doi.org/10.1016/j.foodcont.2016.11.037>

13. To KK, Lu L, Yip CC, Poon RW, Fung AM, Cheng A, Lui DH, Ho DT, Hung IF,

Chan KH, Yuen KY. Additional molecular testing of saliva specimens improves

the detection of respiratory viruses. Emerg Microbes Infect. 2017 Jun 7;6(6):e49.

doi: 10.1038/emi.2017.35

14. Vogels CB, Watkins AE, Harden CA, Brackney D, Shafer J, Wang J, et al.

SalivaDirect: A simplified and flexible platform to enhance SARS-CoV-2 testing

capacity. Med. 2020. <https://doi.org/10.1016/j.medj.2020.12.010>.

15. Watkins AE, Fenichel EP, Weinberger DM, Vogels CBF, Brackney DE,

Casanovas-Massana A, Campbell M, Fournier J, Bermejo S, Datta R, Dela Cruz

CS, Farhadian SF, Iwasaki A, Ko AI, Grubaugh ND, Wyllie AL; Yale IMPACT

- Research Team2. Increased SARS-CoV-2 Testing Capacity with Pooled Saliva Samples. *Emerg Infect Dis.* 2021 Apr;27(4):1184–7. doi: 10.3201/eid2704.204200.
16. World Health Organization (WHO). Weekly epidemiological update on COVID-19 - 1 June 2023. <https://www.who.int/publications/m/item/weekly-epidemiological-update-on-covid-19---1-june-2023> (accessed June 5, 2023).
17. Yamazaki W, Makino R, Nagao K, Mekata H, Tsukamoto K. New Micro-amount of Virion Enrichment Technique (MiVET) to detect influenza A virus in the duck faeces. *Transbound Emerg Dis.* 2019 Jan;66(1):341-348. doi: 10.1111/tbed.13027.
18. Yamazaki W, Matsumura Y, Thongchankaew-Seo U, Yamazaki Y, Nagao M. Development of a point-of-care test to detect SARS-CoV-2 from saliva which combines a simple RNA extraction method with colorimetric reverse transcription loop-mediated isothermal amplification detection. *J Clin Virol.* 2021 Mar;136:104760. doi: 10.1016/j.jcv.2021.104760.
19. Yamazaki Y, Thongchankaew-Seo U, Yamazaki W. Very low likelihood that cultivated oysters are a vehicle for SARS-CoV-2: 2021-2022 seasonal survey at supermarkets in Kyoto, Japan. *Heliyon.* 2022 Oct;8(10):e10864. doi: 10.1016/j.heliyon.2022.e10864.

Table 1. LOD determination of developed and conventional methods using saliva spiked with a 10-fold dilution series of SARS-CoV-2.

Test 1	Virus dilution	Neat	1:10	1:100	1:1000	1:10000
	Number of virus copies spiked per ml of saliva	100,000	10,000	1,000	100	10
	Estimated virus number of copies per PCR reaction tube	28,800	2,880	288	29	3
	Developed method	ND	31.507±0.222	36.078±0.198	No. Ct	No. Ct
	Estimated virus number of copies per PCR reaction tube	160	16	2	0.2	0.02
	Conventional method (PBS)	38.211±0.457	(38.967*)	No. Ct	No. Ct	ND
	Conventional method (SAP)	37.526±0.952	(38.694 ⁺)	No. Ct	No. Ct	ND
Test 2	Virus dilution	Neat	1:10	1:100	1:1000	1:10000
	Number of virus copies spiked per ml of saliva	40,000	4,000	400	40	4
	Estimated virus number of copies per PCR reaction tube	11,520	1,152	115	12	1
	Developed method	ND	36.1627±0.429	38.320±1.083	(39.374*)	No. Ct
	Estimated virus number of copies per PCR reaction tube	64	6	0.6	0.06	0.006
	Conventional method (PBS)	36.104±0.917	(39.666 ⁺)	No. Ct	No. Ct	ND
	Conventional method (SAP)	36.310±1.291	No. Ct	No. Ct	No. Ct	ND
Test 3	Virus dilution	Neat	1:10	1:100	1:1000	
	Number of virus copies spiked per ml of saliva	20,000	2,000	200	20	
	Estimated virus number of copies per PCR reaction tube	5,760	576	58	6	
	Developed method	29.2142±0.180	32.280±0.379	38.089±0.029	No. Ct	
	Estimated virus number of copies per PCR reaction tube	32	3	0.3	0.03	
	Conventional method (PBS)	37.718±0.705	No. Ct	No. Ct	No. Ct	
	Conventional method (SAP)	(38.933 ⁺)	(39.651*)	No. Ct	No. Ct	

LOD, Limit of detection.

448 No. Ct, No threshold cycle values detected using real-time PCR.

449 PBS, Phosphate-buffered saline.

450 RT-qPCR, Reverse transcription quantitative real-time polymerase chain reaction.

451 SAP, Semi-alkaline proteinase.

452 SARS-CoV-2, Severe acute respiratory syndrome coronavirus 2.

453

454 Positive results were determined only if all three Ct values were within 40.000 in triplicate.

455 LODs are shown in bold.

456 The figures in parentheses indicate one or two positive results out of three.

457 *, One positive in triplicate analysis.

458 +, Two positives in triplicate analysis.

459

Table 2. Comparison of LOD between developed and conventional methods by RT-qPCR.

	LOD per ml of saliva	LOD per RT-qPCR reaction tube
Developed	200	58
Conventional (PBS)	20,000	32
Conventional (SAP)	40,000	64

LOD, Limit of detection.

PBS, Phosphate-buffered saline.

RT-qPCR, Reverse transcription quantitative real-time polymerase chain reaction.

SAP, Semi-alkaline proteinase.

Figure legends.

Figure 1. Developed protocol for virion concentration from saliva.

SAP, Semi-alkaline proteinase.

1. Mix 12 ml of saliva with 12 ml of SAP (1:1) and keep for 15 min at room temperature.
2. After centrifugation at 4,000 *g* for 5 min, carefully transfer 18 ml of the supernatant into a new 50-ml tube.
3. Add 15 ml of SAP to 18 ml of the supernatant, mix using a vortex and then keep at room temperature for 15 min.
4. After adding 13.2 ml of PEG-NaCl solution, mix by vortexing.
5. After centrifugation at 8,000 *g* for 20 min, carefully discard the supernatant.
6. Add 100 μ l of PBS and dissolve the invisible precipitates by pipetting and scraping with a 1-ml long tip (10 times each).
7. Place a 1-ml short tip into a 50-ml tube and vortex to completely dissolve the precipitate (supernatant residue after flushing + PBS $\div$ 200 μ l), and then transfer to an RNA extraction tube.

Figure 2. Virus detection from saliva with developed and conventional methods.

SAP, Semi-alkaline proteinase.

Conventional method: Nucleic acid extraction of 200 µl of centrifuged supernatant comprising of saliva 50 µl + PBS or SAP 150uL (4x dilution).

Developed method: Nucleic acid extraction after a simple concentration of 12ml of saliva.

The illustrations are cited from the following sources, all used in compliance with the terms and conditions.

Pipettes and Micropipettes: Irasutoya (irasutoya.com)

Vortex mixer: Kagaku Irasuto (science-illust.com)

50 ml tubes: Kenkyu Net (wdb.com/kenq/illust)

Automatic nucleic acid extractor: Precision System Science Co. Ltd. (pss.co.jp/)

Fig. 1

Developed protocol for virion concentration from saliva


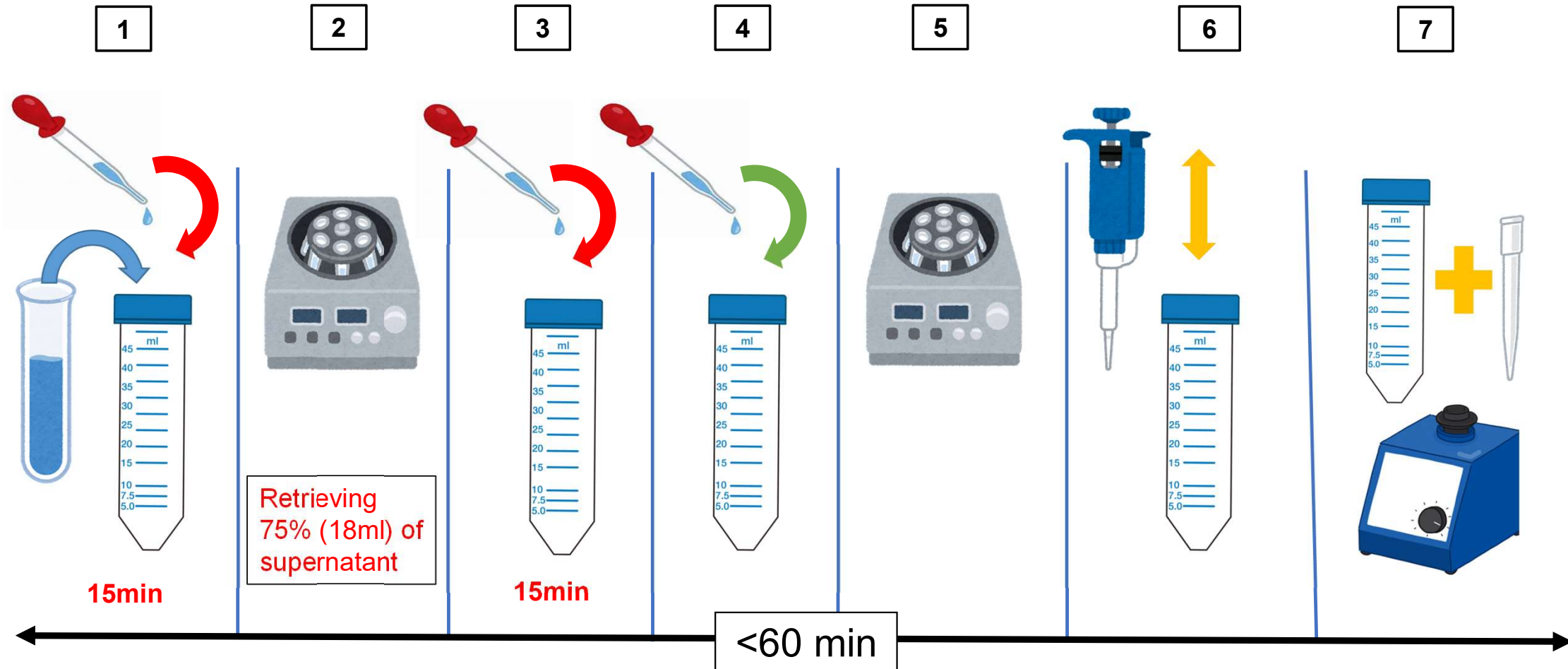


Fig. 2 Virus detection from saliva with developed and conventional methods

Conventional method

50 μ l of saliva

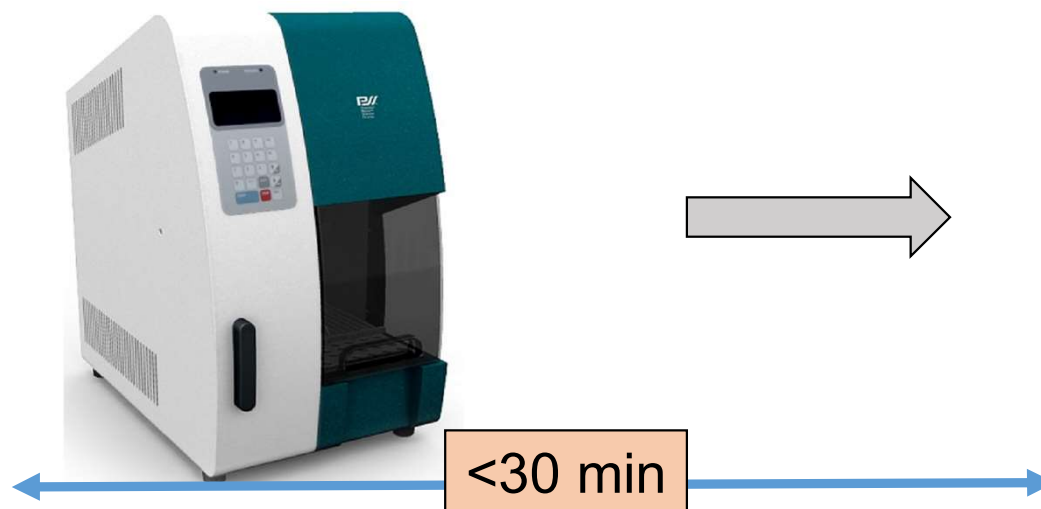

50 μ l of nucleic acid
(1 time concentrated)

180 times $\div$ 7.5 power
difference of 2 (equivalent
to Ct value of 7.5)

50 μ l of nucleic acid
(180 times
concentrated)

Developed method

12ml of saliva

