

SUPPORTING INFORMATION

An integrated lab-on-a-chip device for RNA extraction, amplification and CRISPR-Cas12a-assisted detection for COVID-19 screening in resource-limited settings

Bongkot Ngamsom^a, Alexander Iles^a, Moses Kamita^b, Racheal Kimani^b, Pablo Rodriguez-Mateos^a, Mary Mungai^c, Charlotte E. Dyer^d, Cheryl Walter^d, Jesse Gitaka^{b,1,*}, Nicole Pamme^{a,e,1,*}

^aDepartment of Chemistry and Biochemistry, University of Hull, Hull, HU6 7RX, UK

^bDIRECTORATE OF RESEARCH AND INNOVATION, MOUNT KENYA UNIVERSITY, THIKA, KENYA

^cCentre for Virus Research, Kenya Medical Research Institute, Nairobi, Kenya

^dDepartment of Biomedical Sciences, University of Hull, Hull, HU6 7RX, UK

^eDepartment of Materials and Environmental Chemistry, Stockholm University, SE-106 91, Sweden

¹To whom correspondence should be addressed. Email: nicole.pamme@mmk.su.se; jgitaka@mku.ac.ke.

*Joint senior authors

RNA extraction from viral particles using silica paramagnetic beads and GuHCl solution

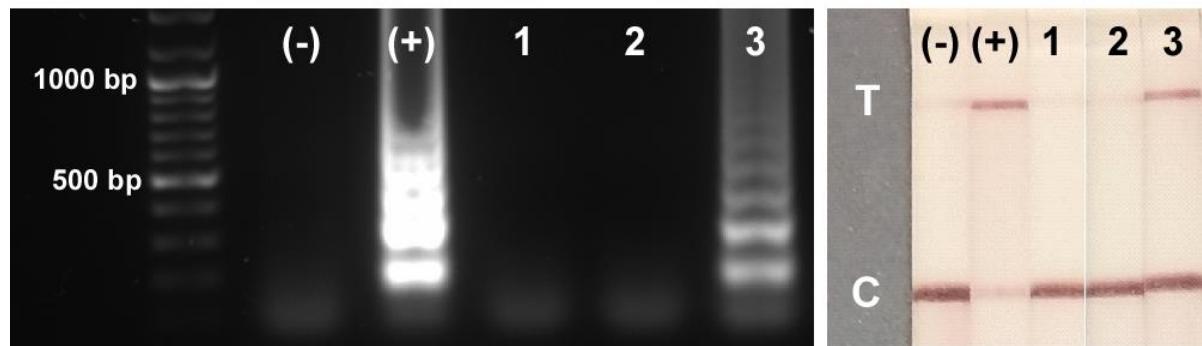


Fig. S1 Gel electrophoresis of tube-based RT-LAMP, and lateral flow test strips from CRISPR-Cas12a detection of RT-LAMP products. Samples: (-) = NTC; (+) = 470 copies of genomic SARS-CoV-2 (free RNA); 1 = 100 copies mL⁻¹ of SARS-CoV-2 verification panel; 2 = MB + SARS-CoV-2 verification panel in nuclease-free water (1000 copies mL⁻¹); 3 = MB + SARS-CoV-2 verification panel in 5 M GuHCl (1,000 copies mL⁻¹). Incubation time for MB and SARS-CoV-2 verification panel = 10 min. RT-LAMP was performed at 64 °C for 35 min using primers targeting N gene. CRISPR-Cas12a assays were carried out at 37 °C for 10 min.

Cigarette filter coupled with IFAST-CRISPR/Cas device for saliva collection and SARS-CoV-2 detection

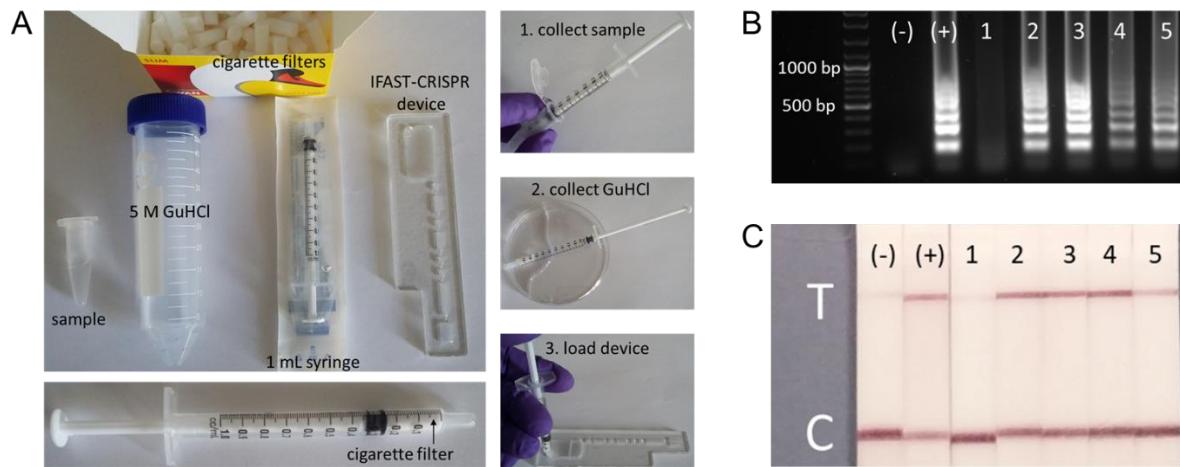


Fig. S2 A) Interfacing of Swan brand cigarette filter for sample collection and IFAST-CRISPR on-chip assays. B) Gel electrophoresis of tube-based RT-LAMP, and C) lateral flow test strips from CRISPR-Cas12a detection of RT-LAMP products. Samples: (-) = NTC; (+) = 470 copies of genomic SARS-CoV-2 (free RNA); 1 = 100 copies mL^{-1} of SARS-CoV-2 verification panel; 2 = MB + SARS-CoV-2 verification panel in nuclease-free water ($1,000$ copies mL^{-1}); 3 = MB + SARS-CoV-2 verification panel in 5 M GuHCl ($1,000$ copies mL^{-1}). Incubation time for MB and SARS-CoV-2 verification panel = 10 min. RT-LAMP was performed at 64°C for 35 min using primers targeting N gene. CRISPR-Cas12a assays were carried out at 37°C for 10 min.

Table S1. Cost breakdown for device and reagents.

	Component	Supplier	Amount per purchase unit		Required amount per device/assay		Number of device/assay per unit	Unit price (USD)	Price per device/assay (USD)
			value	unit	value	unit			
Device materials	PMMA device	Kingston Plastics	1	sheet	0.05	sheet	20	4.70	0.24
	Adhesive PCR plate seals	Thermofisher Scientific	100	sheet	0.5	sheet	200	115.73	0.58
Reagents	Magnesil Beads	Promega	25	mL	0.001	mL	25,000	223.49	0.01
	GuHCl	Promega	100	g	0.478	g	209.36	61.91	0.30
	Milenia HybriDetect 1 Lateral flow strips	TwistDx	100	strip	1	strip	100	296.26	2.96
	LbCas12a	New England Biolabs	70	µL	1	µL	70	87.46	1.25
	gRNA – N gene	Eurofins	122,048	µL	1.25	µL	152,560	603.01	0.004
	gRNA – E gene	Eurofins	122,048	µL	1.25	µL	152,560	262.28	0.0017
	gRNA – RNase P	Eurofins	122,048	µL	1.25	µL	152,560	262.28	0.0017
	LF reporter	Integrated DNA Technologies	2107	µL	0.5	µL	4,214	79.84	0.02
	LAMP primers – N gene FIP	Integrated DNA Technologies	931	µL	0.4	µL	2,328	4.96	0.002
	LAMP primers – N gene BIP	Integrated DNA Technologies	956	µL	0.4	µL	2,390	4.72	0.002
	LAMP primers – N gene F3	Integrated DNA Technologies	10,250	µL	0.05	µL	205,000	2.18	1.06E-05
	LAMP primers – N gene B3	Integrated DNA Technologies	11,000	µL	0.05	µL	220,000	2.66	1.21E-05
	LAMP primers – N gene LF	Integrated DNA Technologies	3,675	µL	0.1	µL	36,750	2.3	6.26E-05
	LAMP primers – N gene LB	Integrated DNA Technologies	4,525	µL	0.1	µL	45,250	2.18	4.82E-05
	LAMP primers – E gene FIP	Integrated DNA Technologies	931	µL	0.4	µL	2,328	4.96	0.002
	LAMP primers – E gene BIP	Integrated DNA Technologies	956	µL	0.4	µL	2,390	4.72	0.002
	LAMP primers – E gene F3	Integrated DNA Technologies	10,250	µL	0.05	µL	205,000	2.18	1.06E-05
	LAMP primers – E gene B3	Integrated DNA Technologies	11,000	µL	0.05	µL	220,000	2.66	1.21E-05

	LAMP primers – E gene LF	Integrated DNA Technologies	3,675	µL	0.1	µL	36,750	2.3	6.26E-05
	LAMP primers – E gene LB	Integrated DNA Technologies	4,525	µL	0.1	µL	45,250	2.18	4.82E-05
	LAMP primers – RNase P - FIP	Integrated DNA Technologies	931	µL	0.4	µL	2,328	4.96	0.002
	LAMP primers – RNase P - BIP	Integrated DNA Technologies	956	µL	0.4	µL	2,390	4.72	0.002
	LAMP primers – RNase P - F3	Integrated DNA Technologies	10,250	µL	0.05	µL	205,000	2.18	1.06E-05
	LAMP primers – RNase P - B3	Integrated DNA Technologies	11,000	µL	0.05	µL	220,000	2.66	1.21E-05
	LAMP primers – RNase P - LF	Integrated DNA Technologies	3,675	µL	0.1	µL	36,750	2.3	6.26E-05
	LAMP primers – RNase P - LB	Integrated DNA Technologies	4,525	µL	0.1	µL	45,250	2.18	4.82E-05
	Bst 2.0 DNA polymerase	New England Biolabs	200	µL	1	µL	200	96.85	0.48
	Warmstart RTx Reverse Transcriptase	New England Biolabs	25	µL	0.5	µL	50	87.46	1.75
	dNTP mix	New England Biolabs	800	µL	3.5	µL	228.57	87.46	0.38
								Total	7.99