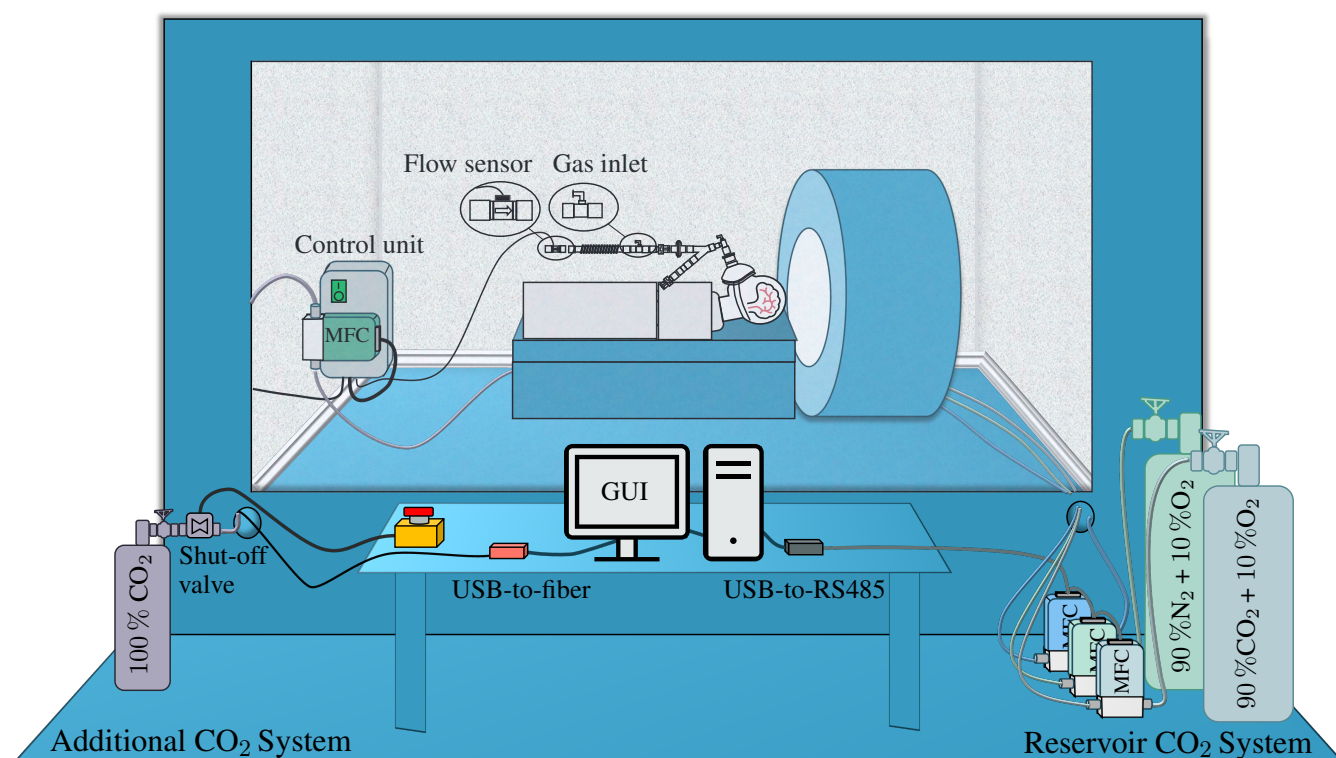


# Supplementary Material

## 1 SYSTEM DESCRIPTIONS AND COMPONENTS

Figure S1 illustrates the configuration of the Additional CO<sub>2</sub> System (on the left) and the Reservoir CO<sub>2</sub> System (on the right) within a magnetic resonance imaging (MRI) environment. In the Additional CO<sub>2</sub> configuration, certain integral components, namely the flow sensor, mass-flow controller, and control unit, are positioned within the confines of the MRI scanner room. In the Reservoir CO<sub>2</sub> configuration, all associated components are situated external to the MRI scanner room.



**Figure S1.** Illustrating the Additional CO<sub>2</sub> and Reservoir CO<sub>2</sub> Systems within an MRI environment. In the Additional CO<sub>2</sub> configuration, the flow sensor, mass-flow controller, and control unit are situated within the confines of the MRI scanner room. This placement choice is required by the necessity to maintain a minimal CO<sub>2</sub> tube length (less than 5 m) between the mass-flow controller and the inlet in the breathing circuit to ensure precise targeting, see primary article figure 2. The mass-flow controller and control unit are affixed to the wall, while the flow sensor is affixed to the breathing circuit, which is positioned on the scanner board. It is worth noting that the flow sensor remains external to the MRI scanner opening, its position adjusted by the expandable tube in the breathing circuit. A computer-based graphical user interface (GUI) and a 100 % CO<sub>2</sub> gas cylinder are both positioned outside the scanner room, connected via tubing and optical cables that traverse through wave guides into the scanner room. Furthermore, the figure illustrates the presence of an emergency shut-off valve, which provides the means for rapid termination of the CO<sub>2</sub> supply. In the Reservoir CO<sub>2</sub> configuration, conversely, all equipment is situated external to the scanner room, with only the gas tubes extending into the room through wave guides to establish connections between the mass-flow controllers and the breathing circuit. The expandable tube in this configuration is extended to a length of 2 m, serving as a reservoir. Notably, the flow sensor is not utilized in the Reservoir CO<sub>2</sub> setup, nor is it affixed to the breathing circuit.

## 1.1 Additional CO<sub>2</sub> System

The Additional CO<sub>2</sub> System is made of various components, which may be categorized into the following main components: CO<sub>2</sub> gas supply, Gas control, Control unit, and the Graphical User Interface (GUI). These components, along with their respective subcomponents, are detailed in table 1.1 and are visually depicted in the left portion of figure S1.

The CO<sub>2</sub> gas supply is made of a gas cylinder with 100 % pure CO<sub>2</sub>, affixed to a shut-off valve, thereby facilitating rapid cessation of gas flow to the mass-flow controller. The Control unit, housing a microcontroller and power supply, controls the setpoint of the mass-flow controller via RS485 protocol and the read the flow sensor via I2C communication to enable a proportional gas flow of CO<sub>2</sub> into the primary respiratory airflow. The required proportional gas flow is calculated at the GUI, given a desired target Fraction of Inspired Carbon Dioxide (FiCO<sub>2</sub>), and subsequently communicated to the Control unit through a fiber optic link.

## 1.2 Reservoir CO<sub>2</sub> System

The Reservoir CO<sub>2</sub> System comprises three fundamental components: a CO<sub>2</sub>/N<sub>2</sub>/O<sub>2</sub> gas supply unit, a gas control system, and a GUI. These constituent elements, along with their respective subcomponents, are detailed in table 1.2 and depicted on the right-hand side of figure S1. The three gas supplies contain 90 % CO<sub>2</sub> + 10 % O<sub>2</sub>, 90 % N<sub>2</sub> + 10 % O<sub>2</sub> and 100 % O<sub>2</sub>, respectively, and are connected to mass-flow controller to enable precise formulation of gas mixtures. The setpoint for each mass-flow controller is calculated at the GUI, given a target concentration and overall gas flow, and communicated over RS485. The outlets of the mass-flow controllers are attached to tubes which and combined using a 4-way connector before being attached to the gas inlet of breathing circuit, see figure 2 in the primary article.

## 1.3 Breathing Circuit

The constituents of the two respiratory circuits employed, depicted in figure 2 of the primary article, are listed in table 1.3.

# 2 BOLD-CVR EXPERIMENT

The final phase of our investigation entailed a test-retest examination of Blood Oxygenation Level Dependent Cerebrovascular Reactivity (BOLD-CVR) with the aim of making a qualitative comparison between the Additional CO<sub>2</sub> and Reservoir CO<sub>2</sub> methods. The experimental configuration is depicted in figure S1.

## 2.1 MRI Protocol

A 3 T Siemens Prisma Magnetic Resonance scanner was employed to monitor alterations in cerebral blood flow resulting from the inhalation of carbon dioxide by the subjects. The MRI protocol encompassed the acquisition of an anatomical T1-weighted MPRAGE scan for image registration. Parameters for this scan were as follows: flip-angle = 8°, TR = 2.3 s, TE = 2.36 ms, TI = 0.9 s, GRAPPA = 3, reference lines = 24, voxel-size = 0.9 mm × 0.87 mm × 0.87 mm, and matrix = 208 × 288 × 288. For each CVR run, utilizing either the Additional CO<sub>2</sub> or Reservoir CO<sub>2</sub> technique, a T2\*-weighted multi-band (MB) gradient-echo planar imaging sequence (EPI), as provided by the Centre for Magnetic Resonance Research (CMRR, University of Minnesota, USA) was collected (Moeller et al., 2010; Feinberg et al., 2010; Xu et al., 2013)). The parameters for this sequence were: flip-angle = 60°, TR = 1 s, TE = 12.60 mm, 29.88 mm, 47.14 mm

**Table S1.** Products and manufacturers of components used in the Additional CO<sub>2</sub> System.

Component		Product, Manufacturer	Comment
CO <sub>2</sub> gas supply	Gas cylinder	S5 (AL), AirLiquide	100 % CO <sub>2</sub> , volume 5 L.
	Pressure regulator	Medireg, Linde Gas	One stage regulator, outlet pressure 4.5 bar.
	High pressure tube	CO <sub>2</sub> hospital tube, Linde Gas	Tube between pressure regulator, shut-off valve and mass-flow controller.
	Quick-connectors	P/N DQSA/DQBA-M-4N-K6-SA, DK-LOK	Quick connectors with checker valves for easy disconnection of the tube between the shut-off valve and mass-flow controller.
Emergency shut-off	3/2 way valve	P/N 126149, Burkert	The 3/2 way valve is connected to outlet of the pressure regulator, and upon breaking the power supply by pressing the stop button, it closes.
	Stop button	P/N XALK178F, Schneider	
Gas control	Mass-flow controller	SFC5500, Sensirion	Calibrated and connected to CO <sub>2</sub> gas source. Max. flow: 10 L min <sup>-1</sup> . Communication via RS485.
	Low pressure tube	Tubclair Al, Tricoflex	Tube between MFC and gas inlet of breathing circuit. ID: 2 mm, OD: 4 mm.
	Flow sensor	SFM3200, Sensirion	Flow range: -100 to 250 L min <sup>-1</sup> . Communication via I2C.
	Quick-connectors	P/N PMC2404/PMCD1004, Colder	Quick connectors with checker valve for easy disconnection of the tube connected to the outlet of the mass-flow controller.
Control unit	Microcontroller	Arduino Beetle, DFRobot	Controlling the mass-flow controller and flow sensor and communication with GUI.
	Software	Arduino program, Gustav Magnusson (author)	In-house developed.
	UART-optic-converter	SPX17510, Sparkfun	To enable fiber communication between Microcontroller and GUI.
	Power supply	RPT-60, Mean Well	Power supply for components in Control unit as well as the mass-flow controller and flow sensor
GUI	Computer	MacBook Pro, Apple	Any personal computer with Python installed can be used as computer.
	Software	Python program, Gustav Magnusson (author)	
	USB-fiber-converter	SPX17508, Sparkfun	The USB converter connects to the computer to enable communication with the Control unit over fiber.
	Fiber optic cable	P/N 164-4255, RS PRO	

and 64.42 mm, MB-factor = 4, GRAPPA = 2, reference lines = 24, voxel-size = 3 mm  $\times$  3 mm  $\times$  3 mm, matrix = 68  $\times$  68  $\times$  44, and #measurement = 345. Additionally, single-band reference (SBRefs) images were saved. Between each CVR run, the subjects were taken out of the scanner and asked to sit up before continuing.

## 2.2 CO<sub>2</sub> Protocol

A simplified CO<sub>2</sub> stimulation protocol was used, comprising a block-stimulus design with an initial baseline period of 1 min, followed by three alternating cycles of 5 % and 0 % CO<sub>2</sub> over 45 s intervals, concluding with an additional minute of baseline recording. To monitor CO<sub>2</sub> and O<sub>2</sub> levels, a Philips Expression MR400 and its IP5 monitor were employed. The end-tidal CO<sub>2</sub> values were made accessible via the RS232 port of IP5 monitor, operating at a data rate of 1 Hz. The recorded end-tidal CO<sub>2</sub> values were subsequently utilized for the computation of CVR maps based on the BOLD signal acquisitions.

## 2.3 Flow Sensor Interference

An unforeseen complication emerged concerning interference of the flow sensor positioned within the MRI scanner. This interference led to disruptions in the communication between the microcontroller and the flow sensor at the onset of scanning. Further investigation revealed that an alteration in the length of the cable connecting the flow sensor to the microcontroller had increased its susceptibility to electromagnetic interference from the MRI scanner. Initial pilot trials employed a shorter 2 m cable without encountering interference issues, but for subsequent scans, a longer 4.7 m cable was employed to prevent people tripping over the cable. Unfortunately, this extended cable length closely approximated a multiple of the wavelength of the MRI scanner's RF-field (2.36 m). By shortening the cable to 4.3 m, most of the disturbances were mitigated, although not eliminated. By carefully securing the cable along the scanner bore, a stable communication link between the microcontroller and the flow sensor was achieved. Two subjects were successfully scanned using this setup. However, when attempting to scan a third subject, the interference issues resurfaced, indicating that the subject's presence also played a role in the interference.

## 2.4 Preprocessing of MRI Data

The MR data was preprocessed using the automated pipeline fMRIPrep 23.1.4 (Esteban et al., 2019). The T1-weighted images were corrected for intensity non-uniformity and skull-stripped. Volume-based spatial normalization to standard space (MNI152NLin2009cAsym, Fonov et al. (2009)) was performed through nonlinear registration with antsRegistration from the ANTs toolbox (Avants et al., 2008). For each of the BOLD runs a reference volume and its skull-stripped version were generated by using the SBRefs. Head-motion parameters with respect to the BOLD reference (six parameters corresponding three rotation and three translation parameters) was estimated before any spatiotemporal filtering using mcflirt of the FSL toolbox (Jenkinson et al., 2002). BOLD runs were then slice-time corrected to and resampled onto their original, native space by applying the transforms to correct for head-motion. A T2\* was estimated from the preprocessed EPI echoes, by voxel-wise fitting the maximal number of echoes with reliable signal in that voxel to a monoexponential signal decay model with nonlinear regression. The calculated T2\* map was then used to optimally combine preprocessed BOLD across echoes following the method described in Posse et al. (1999). The BOLD reference was then registered to the T1w reference using mri\_coreg (FreeSurfer) followed by flirt (FSL) with the boundary-based registration method. The BOLD time-series were finally resampled into standard space (MNI152NLin2009cAsym) utilizing ANTs and the combined BOLD to T1 and T1 to standard space transformation computed earlier.

## 2.5 Generation of CVR Maps

The generation of CVR maps entailed the utilization of brain-masked BOLD data in conjunction with end-tidal CO<sub>2</sub> samples. First, the BOLD data underwent smoothing through the usage of a Gaussian kernel with a full-width-half-maximum of 6 mm. The CO<sub>2</sub> data was then temporally adjusted to maximize its correlation with the global BOLD signal, thereby establishing a solid initial alignment. Subsequently, the cross-correlation between the time series of each voxel and end-tidal CO<sub>2</sub> was computed, and the CO<sub>2</sub> time series was further adjusted on a voxel-by-voxel basis to the maximum correlation. The voxel-specific temporal adjustments were constrained within a range of  $\pm 20$  s, and the Pearson's correlation p-value, lower-thresholded at 0.05, was used to filter out voxels with low correlation.

The next step involved regressing the voxels against the timeshifted, voxel-specific, end-tidal CO<sub>2</sub> data. To account for motion, the previously computed rotation and translation head-motion parameters, including their temporal derivative and power (total of 24 parameters), were included in the regression as confounding timeseries. The resulting end-tidal CO<sub>2</sub> regression coefficient was used to calculate the percentage change in BOLD signal per unit change in end-tidal CO<sub>2</sub>, denoted as  $\% \Delta \text{BOLD} / \text{mmHg}$ , following the methodology elucidated by Liu et al. (2019).

## 3 VENTILATOR TEST CASES

Fourteen ventilator test cases, prescribed by the European standard ISO 80601-2-12:2020, were incorporated for the purpose of assessing the efficacy and functionality of the Additional CO<sub>2</sub> method within the context of mechanical ventilation. These test cases, including test-lung parameters, can be found in table S4.

## REFERENCES

- Avants, B. B., Epstein, C. L., Grossman, M., and Gee, J. C. (2008). Symmetric diffeomorphic image registration with cross-correlation: Evaluating automated labeling of elderly and neurodegenerative brain. *Medical Image Analysis* 12, 26–41. doi:10.1016/j.media.2007.06.004
- Esteban, O., Markiewicz, C. J., Blair, R. W., Moodie, C. A., Isik, A. I., Erramuzpe, A., et al. (2019). fMRIPrep: A robust preprocessing pipeline for functional MRI. *Nature Methods* 16, 111–116. doi:10.1038/s41592-018-0235-4
- Feinberg, D. A., Moeller, S., Smith, S. M., Auerbach, E., Ramanna, S., Glasser, M. F., et al. (2010). Multiplexed Echo Planar Imaging for Sub-Second Whole Brain fMRI and Fast Diffusion Imaging. *PLOS ONE* 5, e15710. doi:10.1371/journal.pone.0015710
- Fonov, V., Evans, A., McKinstry, R., Almli, C., and Collins, D. (2009). Unbiased nonlinear average age-appropriate brain templates from birth to adulthood. *NeuroImage* 47, S102. doi:10.1016/S1053-8119(09)70884-5
- ISO 80601-2-12:2020 (2020). *Medical Electrical Equipment — Part 2-12: Particular Requirements for Basic Safety and Essential Performance of Critical Care Ventilators*. Standard, European Committee for Standardization
- Jenkinson, M., Bannister, P., Brady, M., and Smith, S. (2002). Improved Optimization for the Robust and Accurate Linear Registration and Motion Correction of Brain Images. *NeuroImage* 17, 825–841. doi:10.1006/nimg.2002.1132
- Liu, P., De Vis, J. B., and Lu, H. (2019). Cerebrovascular reactivity (CVR) MRI with CO<sub>2</sub> challenge: A technical review. *NeuroImage* 187, 104–115. doi:10.1016/j.neuroimage.2018.03.047

- Moeller, S., Yacoub, E., Olman, C. A., Auerbach, E., Strupp, J., Harel, N., et al. (2010). Multiband multislice GE-EPI at 7 tesla, with 16-fold acceleration using partial parallel imaging with application to high spatial and temporal whole-brain fMRI. *Magnetic Resonance in Medicine* 63, 1144–1153. doi:10.1002/mrm.22361
- Posse, S., Wiese, S., Gembris, D., Mathiak, K., Kessler, C., Grosse-Ruyken, M. L., et al. (1999). Enhancement of BOLD-contrast sensitivity by single-shot multi-echo functional MR imaging. *Magnetic Resonance in Medicine* 42, 87–97. doi:10.1002/(sici)1522-2594(199907)42:1<87::aid-mrm13>3.0.co;2-o
- Xu, J., Moeller, S., Auerbach, E. J., Strupp, J., Smith, S. M., Feinberg, D. A., et al. (2013). Evaluation of slice accelerations using multiband echo planar imaging at 3 Tesla. *NeuroImage* 83, 10.1016/j.neuroimage.2013.07.055. doi:10.1016/j.neuroimage.2013.07.055

**Table S2.** Products and manufacturers of components used in the Reservoir CO<sub>2</sub> System.

Component		Product, Manufacturer	Comment
CO <sub>2</sub> gas supply	Gas cylinder	K-sized, AirLiquide	90 % CO <sub>2</sub> + 10 % O <sub>2</sub> , volume 50 L.
	Pressure regulator	Multistage, QMT	Two stage regulator, outlet pressure set to 1 bar
	High pressure tube	Super Nobelair, QMT	Tube between pressure regulator and mass-flow controller.
N <sub>2</sub> gas supply	Gas cylinder	K-sized, AirLiquide	90 % N <sub>2</sub> + 10 % O <sub>2</sub> , volume 50 L.
	Pressure regulator	Multistage, QMT	Two stage regulator, outlet pressure set to 1 bar
	High pressure tube	Super Nobelair, QMT	Tube between pressure regulator and mass-flow controller.
O <sub>2</sub> gas supply	Hospital wall outlet	Unknown	100 % O <sub>2</sub>
	Pressure regulator	Gas outlet regulator, QMT	One stage regulator, outlet pressure set to 1 bar.
	High pressure tube	O <sub>2</sub> hospital tube, Linde Gas	Tube between pressure regulator and mass-flow controller.
Gas control	Mass-flow controllers	SLA5850, Brooks	All three controllers are calibrated and connected to respective gas source. Max. flow: 10 L min <sup>-1</sup> . Communication via RS485.
	Low pressure tube	Tygon E3603, Saint-Gobain	Tube between mass-flow controllers and 4-way connector. ID 6.4 mm, OD 9.5 mm.
	4-way connector	P/N X0-4NK, Colder	Connecting the three mass-flow controller tubes to the gas inlet of the breathing circuit
GUI	Computer	MacBook Pro, Apple	Any personal computer with Python installed can be used as computer.
	Software	Python program, Gustav Magnusson (author)	
	USB-RS485-converter	P/N USB-RS485-WE-5000-BT, FTDI Chip	USB converter which enables RS485 communication between computer and mass-flow controllers.



**Table S3.** Products and manufacturers of components used in the Ventilator and Subject Setup breathing circuit.

Component		Product, Manufacturer	Comment
Ventilator Setup	Gas inlet	P/N 1947000 + 2710000, Intersurgical	A luer-connector for the mass-flow controller in Additional CO <sub>2</sub> System to connect to.
	Humidifier	MR290, Fisher & Paykel	An empty humidifier with no water, to mix gas.
	Connector humidifier	P/N 1992000, Intersurgical	Connecting gas inlet connector to humidifier.
	Coaxial ventilator tube	P/N Q712218N-IGS, Envisen	Standard ventilator tube, 2 m long.
	Elbow with sampling port	P/N 2714000, Intersurgical	To sample O <sub>2</sub> /CO <sub>2</sub> gas.
	Test-lung	AccuLung, Fluke Biomedical	Test-lung with variable compliance (10, 20 and 50mL mbar <sup>-1</sup> ) and resistance (5, 20 and 50mbar L <sup>-1</sup> s)
Subject setup	Expandable tube	P/N 1526000, Intersurgical	Expanded in Reservoir CO <sub>2</sub> setup, contracted in Additional CO <sub>2</sub> System.
	Gas inlet	P/N 1947000 + 2710000, Intersurgical	A luer-connector for the mass-flow controllers in Additional CO <sub>2</sub> or Reservoir CO <sub>2</sub> Systems to connect to.
	One-way valve	P/N 1921000, Intersurgical	Two used to control air flow.
	Filter	P/N M1003346, Timik Medical	Particle filter.
	Y-piece	P/N 1901000, Intersurgical	Separates inspiration/expiration side of circuit.
	Outlet adapter	P/N 1967000, Intersurgical	22 mm-female to 22 mm female adapter, to connect outlet of Y-piece to the one-way valve.
	Elbow with sampling port	P/N 2714000, Intersurgical	To sample O <sub>2</sub> /CO <sub>2</sub> gas.
	Face mask adapter	3D-printed, Gustav Magnusson (author)	22 mm-female to 30 mm male adapter.
	Face mask	Mask 7450 V2, Vyaire	Variable sized (XS to L).



**Table S4.** The 14 ventilator test cases used to evaluate the Additional CO<sub>2</sub> method in mechanical ventilation. Note that the *Test case - Number* refers to the tables 201.104/201.105 in the European standard ISO 80601-2-12:2020.

Test case		Test lung parameters		Ventilator settings					
Number	Control mode	Compliance (mL mbar <sup>-1</sup> )	Resistance (mbar L <sup>-1</sup> s)	Frequency (breaths/min)	Insp. time (s)	FiO <sub>2</sub> (%)	PEEP (cmH <sub>2</sub> O)	Insp. pressure (cmH <sub>2</sub> O)	Tidal volume (L)
1	Volume	50	5	20	1	30	5	–	500
2	Volume	50	20	12	1	90	10	–	500
3	Volume	20	5	20	1	90	5	–	500
4	Volume	20	20	20	1	30	10	–	500
5	Volume	20	20	20	1	30	5	–	300
6	Volume	20	50	12	1	90	10	–	300
7	Volume	10	50	20	1	30	10	–	300
1	Pressure	50	5	20	1	30	5	10	–
2	Pressure	50	20	12	1	90	10	15	–
3	Pressure	20	5	20	1	90	5	25	–
4	Pressure	20	20	20	1	30	10	25	–
5	Pressure	20	20	20	1	30	5	15	–
6	Pressure	20	50	12	1	90	10	10	–
7	Pressure	10	50	20	1	90	5	15	–