


Graphical abstract

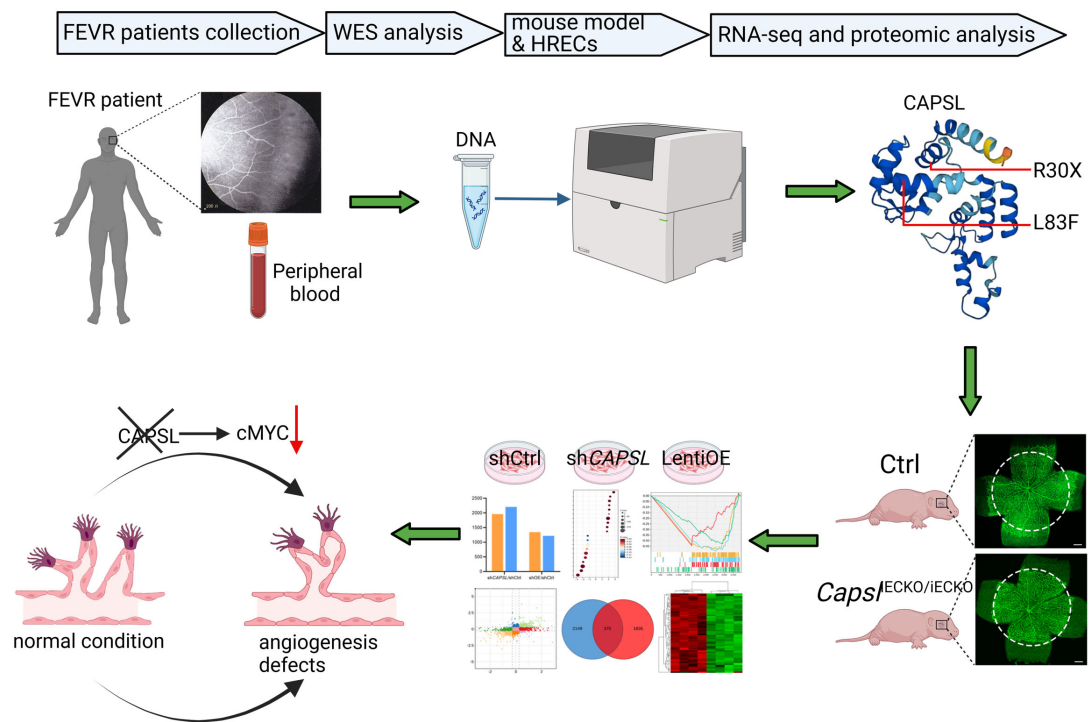

Defective CAPSL function causes impaired retinal angiogenesis through the MYC axis and is associated with familial exudative vitreoretinopathy

Wenjing Liu^{1,2,4#}, Shujin Li^{1,2#}, Mu Yang^{1,2#}, Jie Ma^{1#}, Lu Liu¹, Ping Fei³, Qianchun Xiang¹, Lulin Huang¹, Peiquan Zhao^{3*}, Zhenglin Yang^{1,2*}, Xianjun Zhu^{1,2,4*}

Supplementary data include Figure S1-S2, Table S1 to S3.

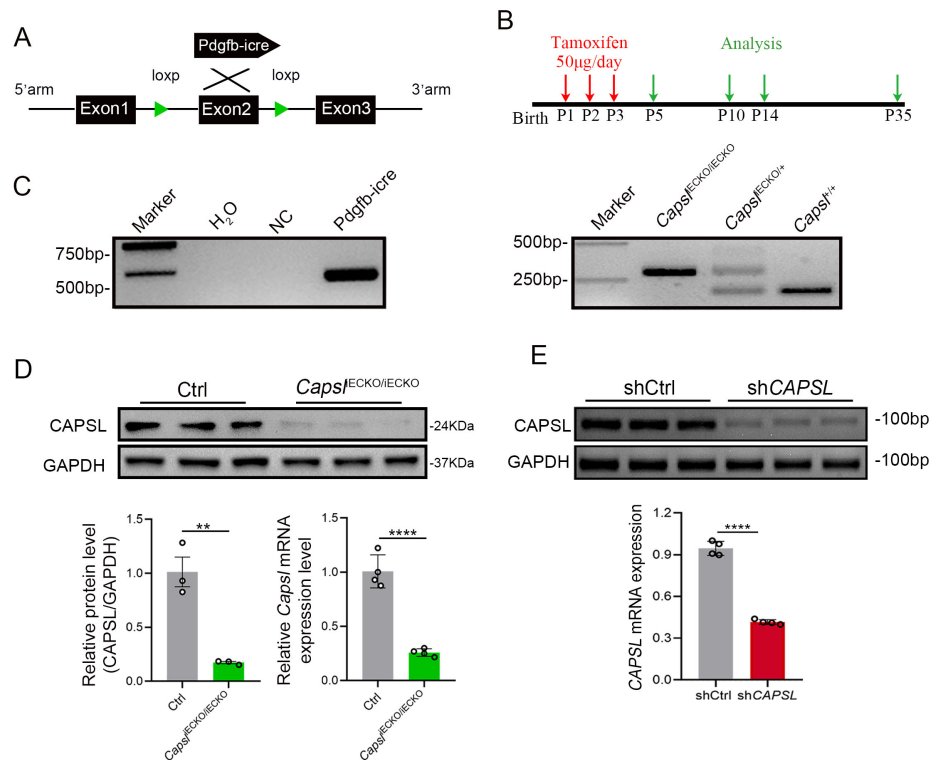


Figure S1. Construction of EC-specific *Capsl* knockout mice model and CAPSL knockdown HREC cell model. (A) Strategy to generate a conditional *Capsl* allele in which exon 2 are knockout specifically in endothelial cells. (B) Diagram depicting the tamoxifen injection and experiment schedule for EC-specific deletion of *Capsl* in retinal vessels from P1 to P3 and their analyses timeline. (C) PCR of genomic DNA from Ctrl and *Capsl*^{IECKO/IECKO} pups. Recombination of the floxed *Capsl* allele occurs only in Tamoxifen-injected animals that are Pdgfrb-iCreER-positive. (D) Immunoblot and quantification of CAPSL knockout efficiency at both protein and mRNA level. Error bars indicate the SD. **P < 0.01, ****P < 0.0001, by Student's t test (n = 3). (E) RT-PCR and RT-qPCR analysis of knockdown efficiency of lentivirus carrying shCAPSL. Error bars indicate the SD. ****P < 0.0001, by Student's t test (n = 3).

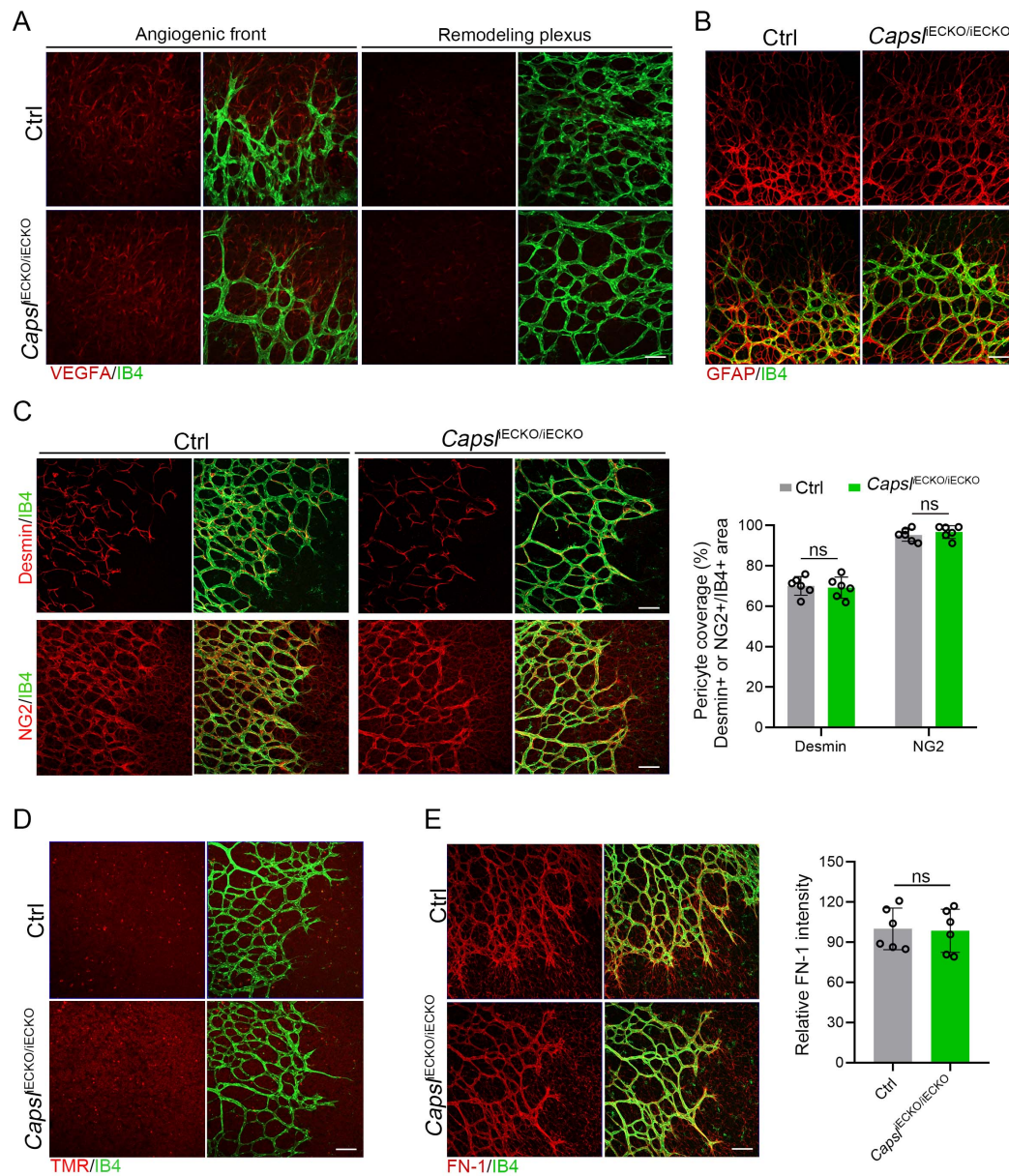


Figure S2. Normal state of other cell types and ECM deposition in *Capsl* knockout mice retina. (A) Co-staining of VEGFA (red) and IB4 (green) at angiogenic front and remodeling plexus in P5 retinas of Ctrl and *Capsl*^{IECKO/IECKO} mice. Scale bar: 50 μ m. (B) Co-staining of Collagen IV (red) and IB4 (green) in P5 retinas of Ctrl and *Capsl*^{IECKO/IECKO} mice. Scale bar: 50 μ m. (C) Maximum intensity projections and quantification of pericyte coverage of the superficial vascular plexus in P5 retinas of Ctrl and *Capsl*^{IECKO/IECKO} mice stained with Desmin/NG2 (red) and IB4 (green). Scale bar: 50 μ m. Error bars indicate the SD. ns: no significance, by Student's t test (n = 6). (D) Representative confocal images of biocytin TMR (red) and IB4 (green) in the retinal flat mounts from P5 Ctrl and *Capsl*^{IECKO/IECKO} mice. Scale bar: 50 μ m. (E) Immunofluorescence staining and quantitative analysis of FN-1 (red) stained ECM deposition in retinas from P5 Ctrl

and *Capsl*^{*iECKO/iECKO*} mice. Scale bar: 50 μ m. Error bars indicate the SD. ns: no significance, by Student's t test (n = 6).

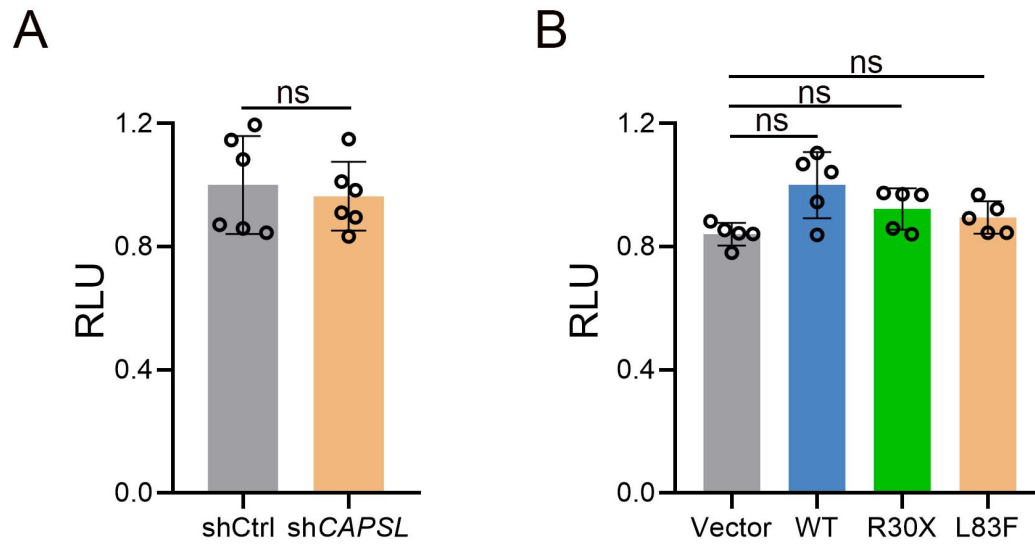


Figure S3. CAPSL does not involve in Norrin/ β -catenin signaling pathway.

(A) Results of luciferase reporter assay of shCtrl-293STF cells and shCAPSL-293STF cells. Error bars indicate the SD. ns: no significance, by Student's t test ($n = 6$). **(B)** Results of luciferase reporter assay of shCAPSL-293STF cells were transfected with plasmids containing CAPSL (WT, R30X, or L83F) or empty vector. Error bars indicate the SD. ns: no significance, by Student's t test ($n = 5$).

Table S1. Pathological prediction of *CAPSL* variants

Software	Variant	Prediction
Mutation Taster	R30X	disease causing: amino acid sequence changed protein features (might be) affected splice site changes
	L83F	disease causing: amino acid sequence changed protein features (might be) affected
Polyphen-2	R30X	-
	L83F	Possibly damaging with a score of 0.777 (sensitivity: 0.85; specificity: 0.92)
PROVEAN	R30X	Score:-13.38, prediction: deleterious
	L83F	Score:-2.51, prediction: deleterious

Table S2. Antibodies used for immunofluorescence studies.

Antibody	Manufacturer	States	Cat. No.	Dilution
Isolectin-B4-Alexa488	Thermo Fisher	USA	I21411	1:200
VEGFA	R & D System	USA	AF-493-NA	1:100
Collagen IV	Millipore Sigma	USA	AB756P	1:100
Fibronectin-1	Sigma	USA	F3648	1:100
NG2	Millipore Sigma	USA	AB5320	1:200
Desmin	Abcam	UK	AB15200	1:100
ERG	Abcam	UK	Ab92513	1:200
GFAP	Cell Signaling Technology	USA	80788	1:200
DAPI	Cell Signaling Technology	USA	4083	1:2000
Donkey anti-Rabbit Alexa594	Invitrogen	USA	A-21207	1:500
Donkey anti-Goat Alexa594	Invitrogen	USA	A-11058	1:500
Donkey anti-Rabbit Alexa647	Invitrogen	USA	A-31573	1:500

Table S3. Primer sets used for RT-qPCR.

Name	Sequence 5'-3'	
<i>Capsl</i> -Mus	Forward	GGACTTGGCAGAGTGTTTCG
	Reverse	TCTGGCTCTGGACATTGGAG
<i>Gapdh</i> -Mus	Forward	TGTGTCCGTCGTGGATCTGA
	Reverse	TTGCTGTTGAAGTCGCAGGAG
<i>CAPSL</i> -Homo	Forward	GCCAAGAAAAAGCTCACCAC
	Reverse	CTTTTCCATGACCACAGCA
<i>MYC</i> -Homo	Forward	AACACACAACGTCTTGGAGC
	Reverse	GCACAAGAGTTCCGTAGCTG
<i>GFP</i> -Homo	Forward	ACGTAAACGGCCACAAGTTC
	Reverse	AAGTCGTGCTGCTTCATGTG
<i>GAPDH</i> -Homo	Forward	CCATGGGTGGAATCATATTGGA

	Reverse	TCAACGGATTTGGTCGTATTGG
--	---------	------------------------

Table S4. Antibodies used for western blot analysis.

Antibody	Manufacturer	States	Cat. No.	Dilution
CAPSL	Proteintech	China	17174-1-AP	1:2000
MCM2	Cell Signaling Technology	USA	3619	1:1000
MCM3	Cell Signaling Technology	USA	4012	1:1000
MCM4	Cell Signaling Technology	USA	3228	1:1000
MCM5	Abcam	UK	AB76023	1:2000
MCM7	Cell Signaling Technology	USA	3735	1:1000
E2F1	Abcam	UK	Ab92513	1:2000
CyclinD1	Cell Signaling Technology	USA	55506	1:2000
CyclinD2	Cell Signaling Technology	USA	3741	1:2000
CyclinD3	Invitrogen	USA	AB_2809647	1:1000
CyclinE1	Cell Signaling Technology	USA	28808	1:1000
CDK4	Cell Signaling Technology	USA	127903	1:1000
PCNA	ABclonal	China	A12427	1:1000
CDC42	Proteintech	China	0155-1-AP	1:1000
RHOA	ABclonal	China	A4855	1:1000
Rac1	ABclonal	China	A7720	1:1000
MYL9	ABclonal	China	A8738	1:1000
P-MYL9-S19	ABclonal	China	AP9055	1:1000
GAPDH	Proteintech	China	1E6D9	1:5000